



# **BTL-I100**

# **BTL-I200**

**USER'S MANUAL**

## CONTENTS

<b>1</b>	<b>USER'S MANUAL .....</b>	<b>3</b>
1.1	Outward appearance of the BTL-1100 and BTL-1200 .....	3
1.2	Configurations of the BTL-1100 and BTL-1200 couches .....	3
1.3	The BTL-1100 couch configuration with folding frame .....	3
1.4	Control units .....	4
1.5	Cleaning and maintenance .....	4
1.6	Safety measures .....	4
1.7	Warranty .....	4
<b>2</b>	<b>TECHNICAL PARAMETERS .....</b>	<b>5</b>
2.1	Technical parameters of the BTL-1100 and BTL-1200 couches .....	5
2.2	Accessories for the BTL-1100 and BTL-1200 couches .....	5
2.3	Classification .....	5
2.4	Procedure of conformity assessment .....	5
2.5	Applicable standards and regulations .....	5

## 1 USER'S MANUAL


### 1.1 OUTWARD APPEARANCE OF THE BTL-1100 AND BTL-1200 COUCHES

The therapy couches BTL-1100 and BTL-1200 consist of:


- firm welded construction
  - fixed height
  - manually adjustable head section tilt
  - the option to divide the main section of the couch into two parts with manually adjustable mid-section elevation
- The couch frame is equipped with height-adjustable plastic leg endings for elimination of surface unevenness. The fixed main section and adjustable head section are made of chipboard with foam rubber, and the both sections are covered with leatherette upholstery. The head section of the couch can be supplied with the so-called "nose cut-away" - removable part for comfortable abdominal position of the patient.

### 1.2 CONFIGURATIONS OF THE BTL-1100 AND BTL-1200 COUCHES

- 2-section couch with manually adjustable head section

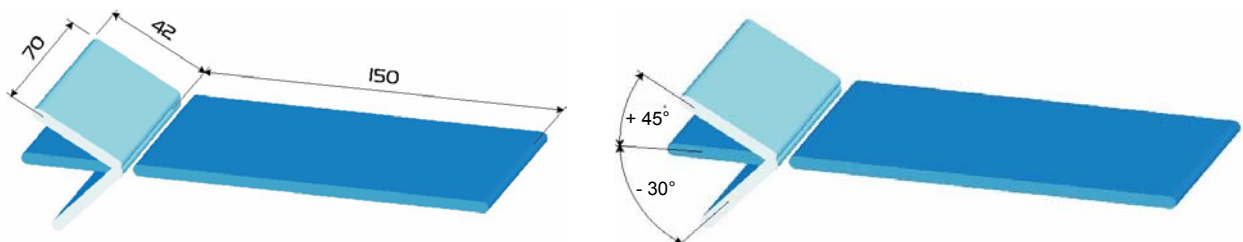

- 3-section couch with manually adjustable head section tilt and mid-section elevation

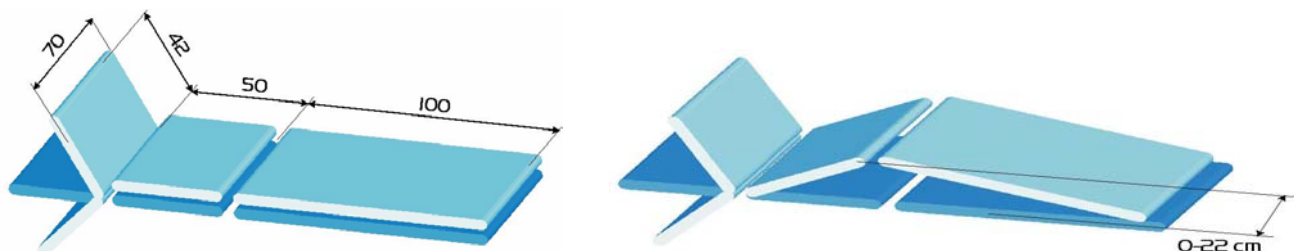

### 1.3 THE BTL-1100 COUCH CONFIGURATION WITH FOLDING FRAME

The BTL-1100 couch with folding frame is the same as the above-given model. The only difference is that during the transport, the spacing bar between couch legs is removed, and thus the legs are folded towards the frame. In this way, the height of the couch is reduced to 20 cm, and it enables its easier transportation and manipulation. When the couch is unpacked, it is necessary to straighten its legs, insert the spacing bar between them, and fix it by the supplied screws.

## **1.4 CONTROL UNITS**

- manual adjustment of the head section tilt (two fixed locking levers)  
To adjust the head section upwards or downwards, push the locking levers towards the head section and release them when the desired angle is reached.
- manual adjustment of the mid-section elevation - 3-section couch configuration (two levers placed in the middle of the couch)  
To adjust the mid-section of the couch, especially downwards, it is necessary to hold this section before releasing the levers.

## **1.5 CLEANING AND MAINTENANCE**

Wash the upholstery with a detergent. For disinfection, use for example 10% Persteril. Do not use acetone-based substances!

## **1.6 SAFETY MEASURES**

When adjusting the head section tilt or mid-section elevation, it is necessary to hold firmly the head section or mid-section before releasing the locking levers to avoid uncontrolled movement of the couch and possible injury of the personnel.

## **1.7 WARRANTY**

The product is warranted against defects in materials and workmanship for a period of 24 months from the date of shipment. In case of any malfunction, contact your BTL distributor. The manufacturer will at its option repair or replace components which prove to be defective during the warranty period provided that the product is used in accordance with the manual. Should the product be returned to the manufacturer for repair or replacement, it must be sent carriage paid.

## 2 TECHNICAL PARAMETERS

### 2.1 TECHNICAL PARAMETERS OF THE BTL-1100 AND BTL-1200 COUCHES

Length: 195 cm	X
Standard width: 70 cm	X
Reduced width: 60 cm	O
Extended width: 80 cm	O
Standard height: 70 cm	X
Height: 60 cm	O
Height: 80 cm	O
Folding frame (only for the BTL-1100)	O
Mid-section elevation: 0 – 22 cm (with the 3-section configuration)	x
Head section adjustment: +45°/ -30°	X
Weight: 25 – 35 kg (according to the model)	X
Load capacity: 170 kg	X

x – standard

o – optional

### 2.2 ACCESSORIES FOR THE BTL-1100 AND BTL-1200 COUCHES

Nose cut-away for the BTL-1100	O
Nose cut-away for the BTL-1200	X
Sanitary paper holder	O
Exercise accessories	O

x – standard

o – optional

### 2.3 CLASSIFICATION

Classification according to statute no. 181/2001 Sb.:	class I
---	---------

### 2.4 PROCEDURE OF CONFORMITY ASSESSMENT

Procedure of conformity assessment:	According to Annex no. XII statute no. 181/2001 Sb. and statute no. 336/2001 Sb.
-------------------------------------	--

### 2.5 APPLICABLE STANDARDS AND REGULATIONS

Standards:	EN ISO 9001, EN 46001, EN 1441
Regulations:	Statute no. 181/2001 Sb. and statute no. 336/2001 Sb.

**Manufacturer:**

Medical Technologies s.r.o.  
Heleny Malirove 11  
169 00 Prague 6  
Czech Republic  
tel/fax: +420 233 355 584

**Representative:**

BTL Industries Limited  
Suite 401 Albany House  
324-326 Regents Street  
London  
W1B 3BL  
United Kingdom

tel: +4420 8371 2267  
fax: +4420 8371 2270  
email: [btl@cobus.net](mailto:btl@cobus.net)

**International Sales:**

BTL s.r.o.  
Heleny Malirove 11  
169 00 Prague 6  
Czech Republic

tel.: +420 220 517 012, 220 514 753  
fax: +420 220 515 520  
e-mail: [sales@btlnet.com](mailto:sales@btlnet.com)  
<http://www.btlnet.com>

© All rights reserved. No part of this manual may be reproduced, saved in a research centre or transferred by any means incl. electronic, mechanic, photographic or other records without previous approval from Medical Technologies s.r.o.

The contents of this document is provided "as is ". Except as required by applicable law, no warranties of any kind, either expressed or implied, are made in relation to the accuracy, reliability or contents of this document. Medical Technologies s.r.o. reserves the right to revise this document or withdraw it at any time without prior notice.